



# ALD

Academy of Laser Dentistry



Experience the Excitement of Lasers  
in Everyday Dentistry  
*A Hands-on Workshop*

Scott D. Benjamin, D.D.S.  
[Technology@DentalAIM.com](mailto:Technology@DentalAIM.com)  
[www.DentalAIM.com](http://www.DentalAIM.com)  
607-563-3927

Buffalo Niagara Dental Meeting  
November 3, 2011  
[www.GLIDM.org](http://www.GLIDM.org)

# **Procedures and Applications Performed with a Soft Tissue Diode Laser**

## **Basic Functions of a Soft Tissue Laser**


- Vaporize (Ablate) Soft Tissue
  - Incise Soft Tissue
  - Excise Soft Tissue
  - Ablate (Erase) Soft Tissue
- Photobiomodulate (Stimulate) Tissue Responses

## **A Partial List of Specific Procedures That Can be Performed with the Assistance of a Soft Tissue Laser**

- Gingivectomy
- Gingivoplasty
- Gingival Troughing
- Periodontal Pocket Debridement Laser Therapy (PDLT)
- Biopsies
- Fibroma Removal
- Implant Uncovering
- Flap Surgery
- Soft Tissue Incisions
- Soft Tissue Excisions
- Destruction of Lesions
- Distal / Proximal Wedge
- Operculectomies
- Excision of Pericoronal Gingiva
- Soft Tissue Crown Lengthening
- Removal of Hyperplastic Tissue
- Exposure of Un-erupted Teeth
- Vestibuloplasty / Frenuloplasty
- Frenectomy / Frenotomy
- Incision and Drainage
- Assisting in Bleaching of Dentition

# Laser Assisted Periodontal Care

## Pocket Debridement Laser Therapy (PDLT)

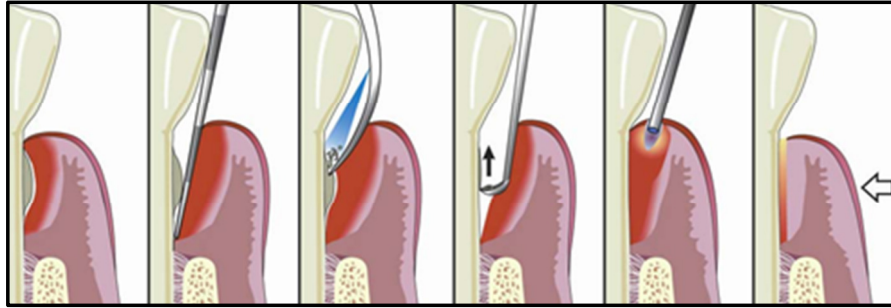

### Overview

When utilizing a laser in periodontal care the historical goals, principles of the causes, treatment, and management of periodontal disease, does not change. The primary goals in the treatment of periodontal disease are to:

- Debride and decontaminate the infected tissue of the periodontal structures
- Decrease the pathogenic bacteria in the periodontal pocket
- Arrest the progressive destruction of periodontal attachment
- Establish a periodontal architecture that can be properly maintained with adequate home care
- Create an environment to facilitate the regeneration of the lost periodontium, whenever possible.

The purpose of laser assisted periodontal care is to restore the periodontal health so that the diseased sulcular and pocket epithelium is eliminated and the pocket (intrasulcular) bacterial activities are minimized. When these goals are accomplished, pocket depths should be minimal and able to be maintained with proper routine homecare, bleeding on probing should be eliminated, and the apical migration of the epithelial attachment halted.

### Treatment Considerations

Comprehensive examination and diagnostic procedures need to be performed to establish an accurate diagnosis of the periodontal status. Prudent clinical judgment is required and the practitioner and the patient need to have realistic expectations for the desired outcomes. Informed consent from the patient on the treatment objectives, prognosis and possibility of tooth loss should be established before any treatment is rendered.

### The Role of the Laser in Laser Assisted Periodontal Care / Periodontal Pocket Debridement Laser Therapy (PDLT)


A soft tissue laser is an adjunctive device primarily used in closed subgingival instrumentation procedures without the displacement (flapping) of the gingiva. These procedures might include root planning and scaling, gingivoplasty, and gingivectomy or a combination of these procedures. The laser's role is not to replace any of the normal procedures or instrumentation, but is used in addition to the normal ultrasonic and hand instrumentation to obtain a better outcome. The laser's primary role is to reduce or eliminate the bacteria along with the diseased, inflamed, or inappropriate soft tissue in the periodontal pocket, and gingiva. This is accomplished by the clinician thoroughly applying photonic energy to the entire soft tissue lining of the periodontal pocket and the extra sulcular diseased epithelium where deemed necessary. The laser also plays a significant role in creating the proper environment for the establishment and organization of a sufficient and stable clot to promote healing. Initiating and maintaining this healing process is imperative in the re-establishment of the periodontal attachment architecture to the root structure of the tooth, thus minimizing pocket depths and to arresting the apical migration of the attachment.

Practitioners need to have a comprehensive understanding of the disease processes, the benefits and limitations of periodontal pocket debridement laser therapy (PDLT) and the present periodontal status of the patient as well as the patient's overall oral and systemic health. The clinician and patient need to understand and remember that periodontal pocket debridement laser therapy (PDLT) is only part of the comprehensive treatment regime, and care that is required to achieve successful outcomes.

# Effects of Light Energy


## Effects of Light on Target Tissue

- Transmission
- Reflection
- Dispersion / Scattering
- Absorption

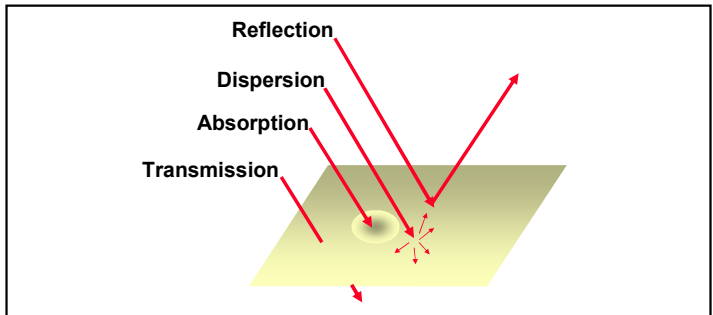

## Effects of Light Energy on Target Subjects

- Photo-thermal Effects
  - Coagulation
  - Vaporization
- Photo-acoustic Effects
  - Disruption
  - Plasma effect
- Light-Induced Fluorescence
  - Caries detection
  - Mucosal evaluation
- Photo-chemical Effects
  - Stimulate chemical reactions
  - Creates chemical bonds
  - Break chemical bonds
- Photobiomodulation
  - Pain relief
  - Wound healing

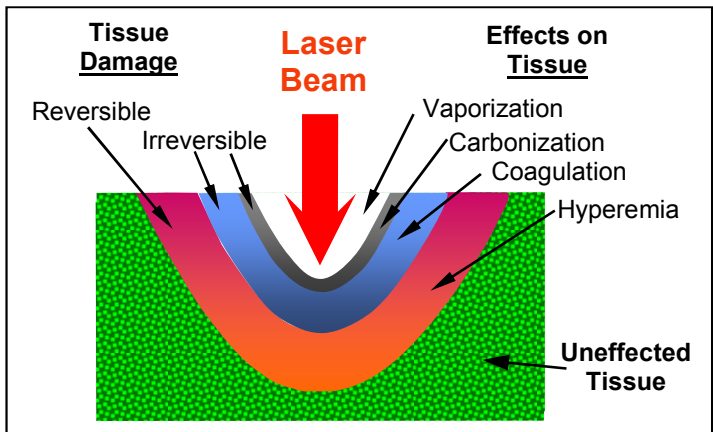

## Thermal Effect of Laser Energy on Tissue

<u>Tissue Temperature (C°)</u>	<u>Observed Effect</u>
37-50°	Hyperthermia
50-60°	Protein Denaturation
> 60°	Coagulation
70-90°	Welding
100-150°	Vaporization
>200°	Carbonization

**There is a linear relationship between the energy of the pulse and the size of the ablation crater.**

Increasing the power lowers the ablation threshold and accelerates the ablation process, thus decreasing thermal side effects.


## Variables Effecting Laser Tissue Interaction

- Wavelength
- Target Composition
  - Chromophores- Substances that absorbs light energy
  - Fluorophores- Substances that emits (produces) light, often when stimulated with light energy
- Interaction Time
  - Temporal Mode
  - Hand Speed
  - Total Interaction Time
- Power
- Energy Transfer Mode
  - Contact vs. Non-Contact
- Spot Size
  - Fiber size (320 micron vs. 200 micron)
- Operator's Knowledge and Experience

## Maximizing the Tissue Interaction Requires:


Matching the proper wavelength with the adequate amount of power with the chromophores present in the tissue.

## Absorption of Laser Energy by Biological Tissue

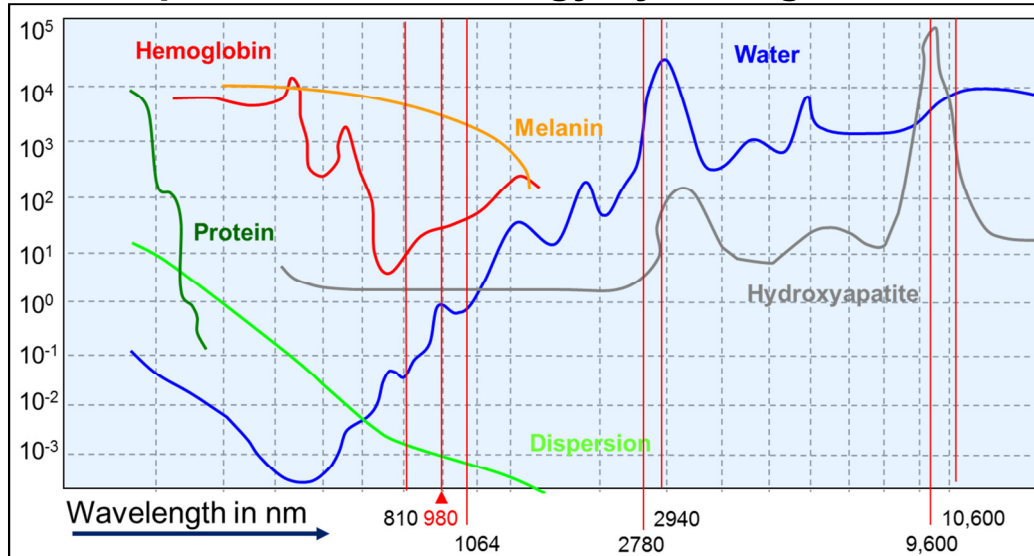

Source: Hale GM, Querry MR, "Optical constants of water in the 200 nm to 200  $\mu$ m wavelength region" Appl. Opt.,12, 555-563

## Transmission & Absorption of Near Infrared (NIR) Light Energy in Water

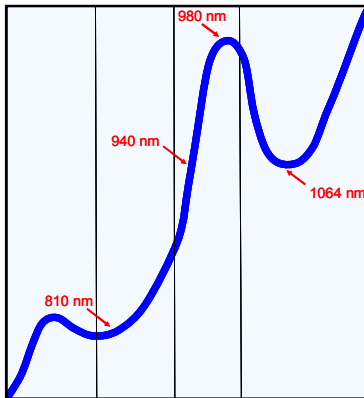

**The Coefficient of Light Transmission in Pure Water is:**


97 % at 810 nm  
85 % at 1064 nm  
65 % at 980 nm

**The Absorption of Light Energy in Water is:**

3% at 810nm  
15% at 1064nm  
35% at 980nm


Cecchetti W, Guazzieri C, Tasca A, Dal Bianco M, Zattoni F, Pagano F  
980 nm Diode Laser and Fiber Optic Resectoscope in Endourological Surgery  
European Biomedical Optics Week, BiOS Europe '96

## Light Absorption of Mucosa for the Near Infrared (NIR) Laser Range

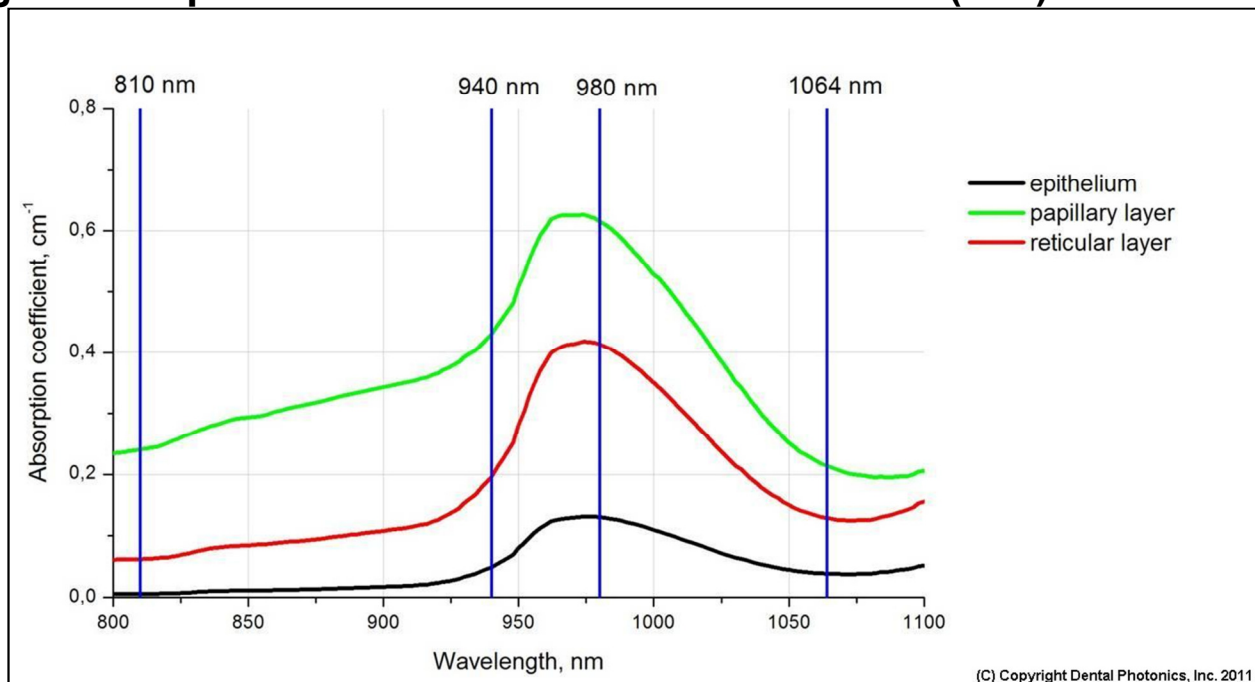

(C) Copyright Dental Photonics, Inc. 2011

## Wavelengths of Light Energy Used in Dentistry on the Electromagnetic Spectrum

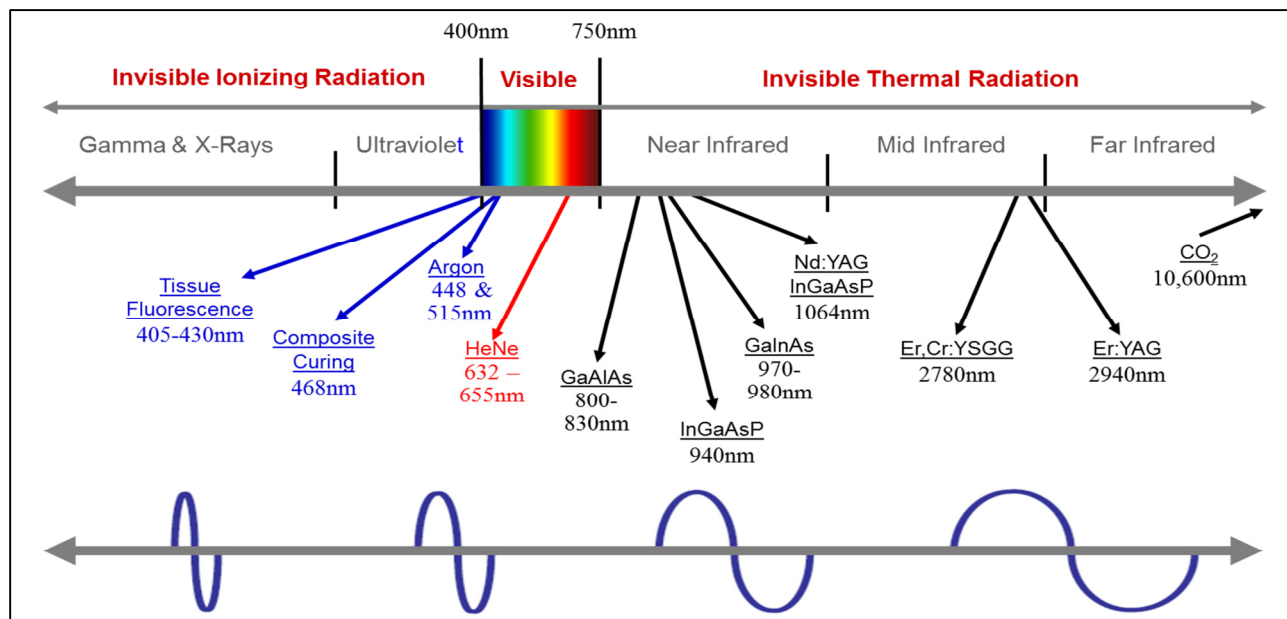

## Active Mediums & Wavelengths of Surgical Dental Lasers

Active Medium	Wavelength(s)
▪ Argon	488-515 nm
▪ Diode	630-1064 nm
• HeNe	630-655 nm
• GaAlAs	805-830 nm
• InGaAsP	940 nm
• GaInAs	970-980 nm
• InGaAsP	1064 nm
▪ Nd:YAG	1064 nm
▪ Erbium,Cr:YSGG	2780 nm
▪ Erbium:YAG	2940 nm
▪ CO <sub>2</sub>	10,600 nm

## Water Content By Percentage (%) in Biological Components

Component / Tissue	Percentage of Water
Mucosa	70%
Skin	70%
Blood	83%
Cartilage	75%
Bone	10-30%
Dentin	12%
Enamel	1-3%



# Temporal Emission Modes

**Temporal Emission Mode** -Managing tissue interaction by controlling the amount of time that laser energy interacts with the tissue, by allowing or not allowing time for the remaining tissue to cool between the pulses of energy being emitted by the laser. The 3 basic temporal emission modes used in dentistry are Continuous Wave (CW), Gated (or Chopped) Pulse, and Free Running Pulse. All the other terms used are variations of the Gated (or Chopped) Pulse mode. Diode lasers have the option of using either a Continuous Wave (CW) or Gated (Chopped) Pulse temporal emission mode.



**Thermal Relaxation Time (TRT)** -The time when the laser energy is not being emitted (or is off).

The purpose of the thermal relaxation time (TRT) is to give the surrounding tissue (the non target tissue) time to cool between the pulses of laser energy thus minimizing the collateral damage and the zones of thermal necrosis and thermal conduction. The longer the TRT the more tissue cooling that occurs.


- **Continuous Wave (CW) Mode** -Does NOT allow for any Thermal Relaxation Time (TRT)
  - CW has the greatest amount of collateral damage and creates the largest zones of thermal necrosis and thermal conduction and therefore has the greatest amount of coagulation.
- **Free Running Pulse Mode** -Provide a long and excellent Thermal Relaxation Time (TRT)
- **Gated (or Chopped) Pulse Mode** - Does allow for Thermal Relaxation Time (TRT)
  - It creates considerably less thermal damage when compared to Continuous Wave (CW).
  - The 980nm diode SIROLaser provides the ability to use water or a liquid for convection cooling and can further reduce the collateral heat spread and damage. When used in a Gated Pulse Mode provides for an even better result by minimizing the zones of thermal necrosis and thermal conduction.

**Duty Cycle (Emission Cycle)** -The temporal duty cycle (often referred to as an emission cycle) is the percentage (%) of time that the laser is Emitting Laser Energy vs. the Thermal Relaxation Time (TRT) within a single pulse. Simply put, it is the percentage of time the laser is on vs. off per pulse cycle.

## Duty Cycle (Emission Cycle) and the Value of Controlling It

An adjustable duty cycle allows the clinician to adjust the duty percentage of time that the energy is being emitted from 1 to 100% of the emission or pulse cycle. This allows for maximum control of the Thermal Relaxation Time (TRT) by extending the TRT as long as desired or even completely eliminating it. This enables the practitioner to accomplish the ideal treatment objective desired and greatly improved outcomes.


## A Comparison of the Thermal Relaxation Time of a 50% and 25% Duty / Emission Cycle

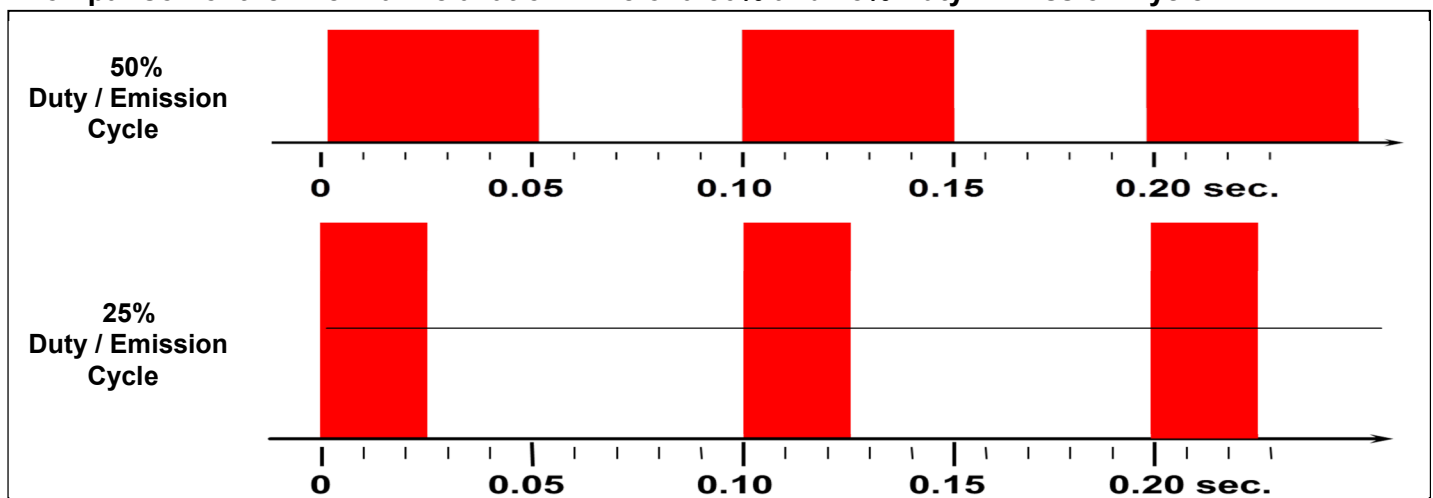


# Temporal Emission Modes

(Continued)


**Continuous Wave Temporal Mode**  
Notice the Continuous Rise in Tissue Temperature  
(Solid Red Line)

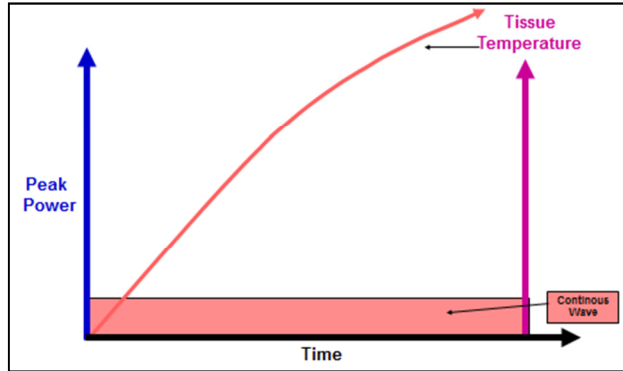


**Gated Pulsed Temporal Mode 50% Duty Cycle**  
Notice the Rise & Fall in Tissue Temperature  
(Zig-Zag Yellow Line)

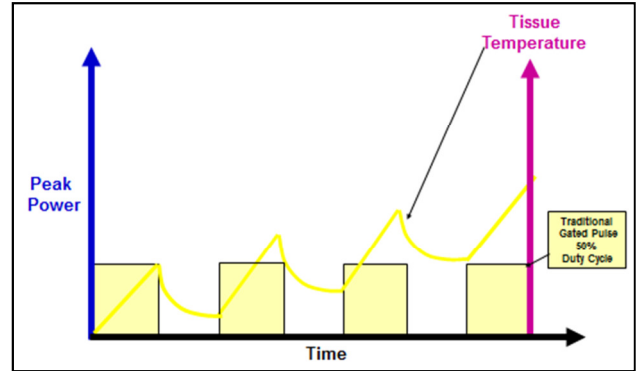

**Gated Pulsed Temporal Mode 12% Duty Cycle**  
Notice the Minimal Rise in Tissue Temperature  
(Zig-Zag Green Line)

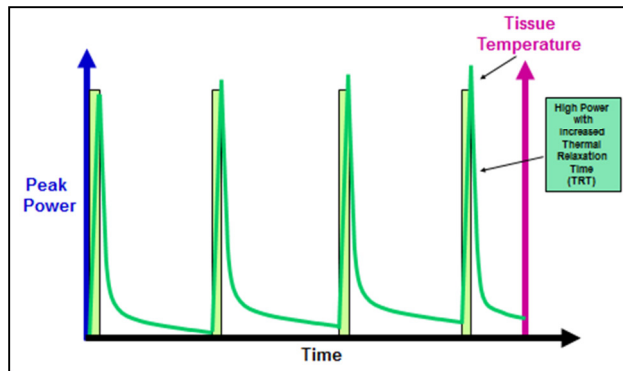

**Comparison of Different Temporal Modes**  
Notice the Difference in Tissue Temperature  
Between Red – Yellow – Green Lines

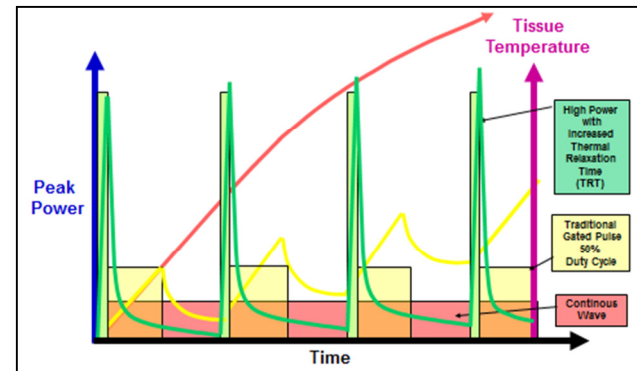

## Comparison of Different Temporal Pulse Modes and Duty / Emission Cycle Effects on Coagulation

Free Running Pulsed  
or Gated Pulsed  
10% Duty Cycle

Gated Pulsed  
25% Duty Cycle

Gated Pulsed  
50% Duty Cycle

Gated Pulsed  
75% Duty Cycle

Continuous  
Wave

No  
Coagulation

Maximum  
Coagulation

Least Amount of  
Collateral Thermal  
Damage

Most Amount of  
Collateral Thermal  
Damage

**Continuous Wave (CW) Mode Will Provide  
Maximum Coagulation and Cause the Most Thermal Damage**



# Laser Device Specifications & Worksheets for Comments & Notes

Use the following pages to record  
your thoughts and comments on your experience with each device

For further information and details on these laser devices you can visit their booths on the exhibit floor  
or contact the companies through their websites and telephone number listed with each device.

For additional information on the use of lasers in dentistry contact the:

## Academy of Laser Dentistry

[www.LaserDentistry.org](http://www.LaserDentistry.org)

Telephone: 954-346-3776

The following are the devices used today in our hands-on workshop:

Laser Model:	<b>Waterlase MD Turbo</b>
Website & Telephone:	<a href="http://www.Biolase.com">www.Biolase.com</a>
Primary Function:	Hard and Soft Tissue Procedures
Wavelength & Medium:	2780 nm Er,Cr:YSGG
Max. Average Power:	8 Watts (Pulse Energy 0-300 mJ)
Nominal Ocular Hazard Distance (NOHD)	2 inches (5 cm)
Temporal Modes:	Free Running Pulsed (10-50 Hz)
Duty Cycle / Pulse Duration:	60 & 700 µsec
Delivery System:	Optical fiber
Tip Design:	Multi use and single use tips
Control Mechanism:	Wireless Foot Control
Powered By:	110v AC
Training / Education:	Several options available
Warranty:	1 Year
Purchased Through:	Henry Schein, or Benco, or Direct from Manufacturer

Comments:

Mfg: **Biolase Technology, Inc.**,  
800-645-6594



Laser Model: **ezLase**  
Website & Telephone: [www.Biolase.com](http://www.Biolase.com)  
Primary Function: All Soft Tissue Procedures,  
Wavelength & Medium: 940 nm Diode  
Max. Peak Power: 7 Watts  
Nominal Ocular Hazard Distance (NOHD) 38 feet 8 inches (11.8 meters)  
Temporal Modes: Continuous Wave and Pulsed  
Duty Cycle: Variable  
Delivery System: Optical fiber  
Tip Design: Single Use Tip  
Control Mechanism: Wireless Foot Control  
Powered By: 110v AC  
Training / Education: DVD and Available for purchase  
Warranty: 2 Years  
Purchased Through: Henry Schein, or Benco, or Direct from Manufacturer

Comments:

Mfg: **Biolase Technology, Inc.,**  
800-645-6594



Laser Model: **iLase**  
Website & Telephone: [www.Biolase.com](http://www.Biolase.com)  
Primary Function: All Soft Tissue Procedures  
Wavelength & Medium: 940 nm Diode  
Max. Peak Power: 5 Watts  
Nominal Ocular Hazard Distance (NOHD) 8 feet 7 inches (2.61 meters)  
Temporal Modes: Continuous Wave and Pulsed  
Duty Cycle: Variable  
Delivery System: Optical Fiber  
Tip Design: Single Use Tip  
Control Mechanism: Finger switch  
Powered By: Rechargeable Li-Ion Battery  
Training / Education: DVD and Available for purchase  
Warranty: 2 Years  
Purchased Through: Henry Schein, or Benco, or Direct from Manufacturer

Comments:

Mfg: **Biolase Technology, Inc.,**  
800-645-6594



Laser Model: **SIROLaser Advance**  
 Website & Telephone: [www.SIROLaserTraining.com](http://www.SIROLaserTraining.com)  
 Primary Function: All Soft Tissue Procedures  
 Wavelength & Medium: 970 nm Diode  
 Max Peak Power: 14 Watts  
 Nominal Ocular Hazard Distance (NOHD) 5 Feet (1.5 meters)  
 Temporal Modes: Continuous Wave & Pulsed  
 Duty Cycle: Variable from 1-100%  
 Delivery System: Optical fiber (320 micron & 200 micron)  
 Tip Design: Single Use Canula & Reusable Autoclavable Fiber  
 Control Mechanism: Finger Switch (Optional Wireless Foot Pedal)  
 Powered By: Rechargeable Battery Powered (4 Hours) or 110v AC  
 Training / Education: 8 Credit Hour Hands-On Device Specific Training Included  
 Warranty: 2 year  
 Purchased From: Major Dental Dealers

Comments:

Mfg: **Sirona Dental Systems**  
 800-659-5977



Laser Model: **DioDent Micro 980**  
 Website & Telephone: [www.ConBio.com](http://www.ConBio.com)  
 Primary Function: All Soft Tissue Procedures  
 Wavelength & Medium: 980 nm Diode  
 Max Peak Power: 3.5 Watts  
 Nominal Ocular Hazard Distance (NOHD) 13 Feet  
 Temporal Modes: Continuous Wave & Pulsed  
 Duty Cycle: 50% Fixed Duty Cycle  
 Delivery System: Optical fiber 300 or 400 microns  
 Tip Design: Disposable Canulas  
 Control Mechanism: Tethered foot control  
 Powered By: 110 vAC  
 Training / Education: Online and free full day off site training  
 Warranty: 2 years  
 Purchased From: Direct from manufacturer

Comments:

Mfg: **Hoya ConBio, Inc.,**  
 800-532-1064



Laser Model: **SL3 – Soft Tissue Laser**  
 Website & Telephone: [www.discusdental.com/lasers](http://www.discusdental.com/lasers)  
 Primary Function: All Soft Tissue Procedures  
 Wavelength & Medium: 808 nm Diode  
 Max Peak Power: 3 Watts  
 Nominal Ocular Hazard Distance (NOHD): 4 feet 8 inches (1.43 Meters)  
 Temporal Modes: Continuous Wave & Pulsed at 10Hz  
 Duty / Emission Cycle: Fixed 50% Duty Cycle (50 msec pulse in pulsed mode)  
 Delivery System: Optical fiber (400 micron)  
 Tip Design: Single Use 400 Micron Fiber Tip  
 Control Mechanism: Wireless Foot Pedal (2.4 GHz RF)  
 Powered By: 9 VDC +/- 5%, 3.5A or Rechargeable Li-Ion Battery  
 Training / Education: Full day 8-CE course through Advanced Laser Training, Inc  
 Warranty: 2 Years  
 Purchased From: Direct from Manufacturer

Mfg: **Discus Dental**

800-217-8822



Laser Model: **NV Microlaser**  
 Website & Telephone: [www.discusdental.com/lasers](http://www.discusdental.com/lasers)  
 Primary Function: All Soft Tissue Procedures  
 Wavelength & Medium: 808 nm Diode  
 Max Peak Power: 2 Watts  
 Nominal Ocular Hazard Distance (NOHD): 3 feet (90.2cm)  
 Temporal Modes: Continuous Wave & Pulsed at 10Hz  
 Duty / Emission Cycle: Fixed 50 % Duty Cycle (50 msec pulse in pulsed mode)  
 Delivery System: Optical fiber (400 micron)  
 Tip Design: Single Use 400 Micron Fiber Tip  
 Control Mechanism: Wireless Foot Pedal (2.4 GHz RF)  
 Powered By: Rechargeable Lithium Ion Battery (Approximately 10 short procedures on full charge)  
 Training / Education: Full day 8-CE course through Advanced Laser Training, Inc  
 Warranty: 1 Year  
 Purchased From: Direct from Manufacturer

Mfg: **Discus Dental**

800-217-8822



Comments:

**Laser Model:** **Odyssey Navigator**  
**Website & Telephone:** [www.IvoclarVivadent.com/odyssey](http://www.IvoclarVivadent.com/odyssey)  
**Primary Function:** All Soft Tissue Procedures  
**Wavelength & Medium:** 810 nm Diode  
**Max Peak Power:** 3 Watts  
**Nominal Ocular Hazard Distance (NOHD)** 13 Feet (154.55 inches)  
**Temporal Modes:** Continuous Wave & Pulsed  
**Duty Cycle:** 50% Fixed Duty Cycle  
**Delivery System:** Optical fiber (400 micron)  
**Tip Design:** Single Use 400 Micron Fiber Tip in 3.5mm or 6mm lengths  
**Control Mechanism:** Wireless Foot Pedal  
**Powered By:** Rechargeable Battery Powered (45 minutes) or 110v AC  
**Training / Education:** Full-day training session at an Ivoclar Vivadent training center  
**Warranty:** 2 years for laser, 1 year for battery, 6 months for fiber optic handpiece  
**Purchased From:** Major Dental Dealers

Mfg: **Ivoclar Vivadent**

800-533-6825



Comments:

**Laser Model:** **Odyssey 2.4G Diode Laser**  
**Website & Telephone:** [www.IvoclarVivadent.com/odyssey](http://www.IvoclarVivadent.com/odyssey)  
**Primary Function:** All Soft Tissue Procedures  
**Wavelength & Medium:** 810 nm Diode  
**Max Peak Power:** 5 Watts  
**Nominal Ocular Hazard Distance (NOHD)** 10 Feet (119.72 inches)  
**Temporal Modes:** Continuous Wave & Pulsed  
**Duty Cycle:** 50% Fixed Duty Cycle  
**Delivery System:** Optical fiber (400 micron)  
**Tip Design:** Disposable Canulas  
**Control Mechanism:** Wireless Foot Pedal  
**Powered By:** 110v AC  
**Training / Education:** Full-day training session at an Ivoclar Vivadent training center  
**Warranty:** 2 years  
**Purchased From:** Major Dental Dealers

Mfg: **Ivoclar Vivadent**

800-533-6825



Comments: